



# DOUBLE DISTILLATION WATER STILL



By **YORCO**  
Scientific Udyog Pvt. Ltd.

**STERI** has been at the forefront of manufacturing Sterilization and Infection control products for over **60 years**. With wide business base, complete in-house design and manufacturing facility, we have developed a reputation for our quality and reliability. At Steri we provide excellent products and exceptional customer service to our clients. Our experience with over **50,000** installations and strive for excellence has resulted in **STERI** being a pioneer in the field.

### Rust-Resistant Construction

High-Quality Stainless Steel is used to avoid corrosion and contamination

### Air Tight Sealing

Joint less steam resistant silicon gasket for effective sealing and safety

### Green Energy Design

A special system is designed to reuse heated cooling water that considerably reduces the energy consumption resulting in saving electricity



### Unique Condenser System

The condenser design ensures that the droplets of distilled water remain in contact with the cooling coil for the longest possible time ensuring maximum energy transfer.

### High Purity Water

Water Distillation equipment provides pyrogen free distillate

### Narrow and Compact

The slim design allows the machine to fit in a small space for maximum area utilization

### Steri Eco Aqua-Save Technology

Steri water distillation system is equipped with water-saving technology to save every drop of water being utilized

### Advanced Sensor System

Steri Water Distillation comes equipped with state of the art smart electronic water sensor for accurate detection and heater protection by controlling the water

## Capacity Available

## Accessories

	VOLUME	LOAD
1	5 Litre/ hour	8 KW
2	10 Litre/ hour	16 KW
3	20 Litre/ hour	28 KW
4	40Litre/ hour	56 KW

Voltage 220V, 50 HZ  
Phase AC

### 1. INCLUDED WITH THE PRODUCT

- Supplied with cord and plug to work on 220/230 volts A.C
- Sturdy Steel Stand (Powder Coated)

### 2. AVAILABLE FOR PURCHASE

- Heater element
- Joint less silicon gasket





## Other Features

All units are made as per IS : 3830 (ISI), CE, ISO: 9001, and ISO:13485 standards

Each unit is fitted with ISI marked immersion heating elements

The instrument is equipped with automatic low water protection system and an easy to use electrical control box

The stand is powder coated for rust free operation and lustrous appearance.

Fuse and current leakage breaker is provided for safety in case of voltage fluctuation or short circuit.

All the electrical and plumbing fixtures are assembled in our factory for quick and easy installation on site

Double Distillation Water Still comes with one year manufacturer's warranty

### Customization of Products :

- Double-walled SS construction with insulation available upon request
- Construction in Stainless Steel 316 available on request
- Fully Automatic control system
- In-built water softner system


## Why Us ?



Afordable



Low Maintenance

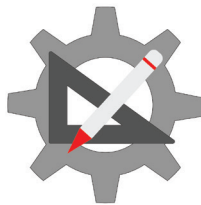

Customizable



Water Saving



Hydraulically Tested at 2.5x Pressure



Multiple Quality Check




1 Year Warranty



Pre - Sale Service



Post - Sale Service

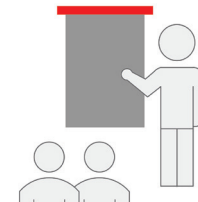

Free Onsite Training



Fast Shipping



Free Installation

## Our Esteemed Clients



## Get In Touch

Marketed by



📍 : D-149 Preet Vihar  
New Delhi 110092.

✉ : [ysupltd@gmail.com](mailto:ysupltd@gmail.com)

🌐 : [yorcoscientific.com](http://yorcoscientific.com)

☎ : 8595293382

Manufactured by



YORCO Scientific Udyog Pvt.Ltd

📍 : 13 Loni road, site 2,  
Mohan nagar, Ghaziabad,  
Uttar Pradesh 201007.

✉ : [ysupltd@gmail.com](mailto:ysupltd@gmail.com)

🌐 : [yorcoscientific.com](http://yorcoscientific.com)

☎ : 9350000607